

FOOD SAMPLING - SOLID SUBSTANCES AND BULK GOODS



Do you want to know what's on your plate?

So do your customers!

Food scandals show us almost daily that food sampling is not considered sufficiently in many enterprises, particularly when inspecting incoming goods. Samples must be taken regularly, because products can only be as good as the raw materials they are made of.

You can carry out very precise material analyzes using the most modern techniques. However, if a mistake is made in sampling then even these modern techniques are useless because the initial sampling error is carried through all the steps of the analysis and incorrect results will be obtained. Thus is the conclusion:

Take samples – but do it right!

Zone samplers

All-layer sampler

For all-layer samples through all the layers of the sample.

Multi sampler

With several closed chambers, for multi point sampling from several predetermined sampling depths.

Uno sampler

Made of stainless steel AISI 316L (1.4404)/PTFE with a closed sample chamber for point sampling from a predetermined sampling depth.





Zone samplers,	Zone samplers, Ø 25 mm (0.98 in.)							
Model	Material	Total length cm (in.)	Insertion depth mm (in.)	No. of openings	Chamber volu- me ml (oz.)	Item no.		
All-layer sampler	AISI 316L	55 (21.65)	430 (16.93)	3	70 (2.37)	5316-2055	N4	
All-layer sampler	AISI 316L	85 (33.46)	710 (27.95)	3	130 (4.40)	5316-2085	N4	
All-layer sampler	AISI 316L	150 (59.06)	1355 (53.35)	5	260 (8.79)	5316-2150	N4	
Multi sampler	AISI 316L	55 (21.65)	430 (16.93)	3	14 (0.47)	5316-3055	N4	
Multi sampler	AISI 316L	85 (33.46)	710 (27.95)	3	17 (0.57)	5316-3085	N4	
Multi sampler	AISI 316L	150 (59.06)	1355 (53.35)	5	17 (0.57)	5316-3150	N4	
Multi sampler	AISI 316L/PTFE	55 (21.65)	430 (16.93)	3	14 (0.47)	5316-0055	N4	
Multi sampler	AISI 316L/PTFE	85 (33.46)	710 (27.95)	3	17 (0.57)	5316-0085	N4	
Multi sampler	AISI 316L/PTFE	150 (59.06)	1355 (53.35)	5	17 (0.57)	5316-0150	N4	
Uno sampler	AISI 316L/PTFE	55 (21.65)	430 (16.93)	1	17 (0.57)	5316-1055	N4	
Uno sampler	AISI 316L/PTFE	85 (33.46)	710 (27.95)	1	17 (0.57)	5316-1085	N4	
Uno sampler	AISI 316L/PTFE	150 (59.06)	1355 (53.35)	1	17 (0.57)	5316-1150	N4	

close-it food & pharma

close-it food & pharma reliably prevents products from leaking out after sampling.

According to FDA guidelines





			_			
close-it food & pha	arma, p	rinted				
W x H mm (in.)	Color	Pieces per roll	Item no.			
95 x 95 (3.74 x 3.74)	red	500	5303-2017	N4		
95 x 95 (3.74 x 3.74)	yellow	500	5303-2020	N4		
95 x 95 (3.74 x 3.74)	blue	500	5303-2021	N4		
95 x 95 (3.74 x 3.74)	green	500	5303-2022	N4		
95 x 95 (3.74 x 3.74)	black	500	5303-2023	N4		
150 x 150 (5.91 x 5.91)	red	250	5303-3017	N4		
150 x 150 (5.91 x 5.91)	yellow	250	5303-3020	N4		
150 x 150 (5.91 x 5.91)	blue	250	5303-3021	N4		
150 x 150 (5.91 x 5.91)	green	250	5303-3022	N4		
150 x 150 (5.91 x 5.91)	black	250	5303-3023	N4		
close-it food & pha	arma, u	inprinted	l			
95 x 95 (3.74 x 3.74)	white	500	5303-2018	N4		
150 x 150 (5.91 x 5.91)	white	250	5303-3018	N4		
close-it ^{tape} food & pharma, continuous roller						
Width mm (in.)	ength m	(ft.)	Item no.			
95 (3.74)	50 (164.0	4)	5303-2019	N4		
150 (5.91)	50 (164.0	4)	5303-3019	N4		



Tubus



Sampler for pourable bulk goods up to approx. 1 cm (0.39 in.) in diameter. Corresponds to ISTA standard for seed lengths of 10 mm (0.39 in.) (wheat, oats) and 5 mm (0.20 in.) (flax, vetches).

- ▶ Also for coarse-grained bulk goods
- ▶ Made of stainless steel AISI 316L (1.4404)
- ▶ Special geometry of the tip

Length cm (in.)	Dia. mm (in.)	Item no.	
40 (15.75)	40 (1.57)	5382-0040	N4
80 (31.50)	40 (1.57)	5382-0080	N4

Milky | acc. to EN ISO 707



Universal powder and granules sampler made of stainless steel AISI 316L (1.4404). Developed initially for milk powder samples, conforms to standards of DIN and IDF (International Dairy Federation).

Volume ml (oz.)	Item no.	
150 (5.07)	5321-0300 N	4

QuickPicker



The core sampler is particularly suitable in the event that bulk goods are taken directly from sacks or bags and have to be filled into the original sampling bottle directly. QuickPicker is delivered complete with two 250 ml (8.45 oz.) PE sample bottles and with a cleaning brush. Corresponds to ISTA standard for seed lengths of 10 mm (0.39 in.) (wheat, oats) and 5 mm (0.20 in.) (flax, vetches).

Material	Insertion depth mm (in.)			Item no.	
AISI 316L	300 (11.81)	50 (19.69)	75 (2.54)	5322-1300	N4

MultiDispo








Disposable zone samplers for taking samples from multiple layers of the sample material. Suitable for free-flowing powders and small-grain granulates.

- ▶ HDPE
- Clean room manufactured, individually packaged

Length mm (in.)	Volume ml (oz.)	Model	Pieces per pack	Item no.	
500 (19.69)	100 (3.38)	standard	20	5393-4440	N4
500 (19.69)	100 (3.38)	sterilized	20	5393-4441	N4
1000 (39.37)	200 (6.76)	standard	20	5393-4450	N4
1000 (39.37)	200 (6.76)	sterilized	20	5393-4451	N4

PowderDispo



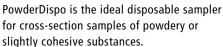

- ▶ HDPF
- Clean room manufactured, individually packaged

Length mm (in.)	Volume ml (oz.)	Model	Pieces per pack	Item no.	
300 (11.81)	60 (2.03)	standard	20	5393-5510	N4
300 (11.81)	60 (2.03)	sterilized	20	5393-5511	N4
500 (19.69)	100 (3.38)	standard	20	5393-5522	N4
500 (19.69)	100 (3.38)	sterilized	20	5393-5521	N4

DispoLance





Suitable for free-flowing powders and small-grain granulates. The sample flows through the open inner tube directly into the container held underneath.

- ▶ HDPF
- Clean room manufactured, individually packaged

Length mm (in.)	Model	Pieces per pack	Artikel-Nr.	
300 (11.81)	standard	20	5393-5560	N4
300 (11.81)	sterilized	20	5393-5561	N4
500 (19.69)	standard	20	5393-5570	N4
500 (19.69)	sterilized	20	5393-5571	N4

FoodScoop

Very stable hand scoop for various applications, particularly suited for the food industry. Seamless connection of scoop and hollow handle, therefore without crevices and edges. Especially smooth, highly polished surface.

- ▶ Especially sturdy model
- > Stainless steel AISI 304 (1.4301)
- ▶ With handy suspension eyelet



Capacity ml (oz.)	Total length mm (in.)	L x W mm (in.)	Item no.	
100 (3.38)	215 (8.46)	120 x 60 (4.72 x 2.36)	5371-0100	N3
200 (6.76)	230 (9.06)	140 x 75 (5.51 x 2.95)	5371-0200	N3
250 (8.45)	245 (9.65)	155 x 85 (6.10 x 3.35)	5371-0250	N3
350 (11.83)	290 (11.42)	180 x 100 (7.09 x 3.94)	5371-0350	N3
500 (16.91)	350 (13.78)	220 x 130 (8.66 x 5.12)	5371-0500	N3



SiloDrill | for sampling from silo vehicles

If bulk goods are transported over long distances (in trucks, railway goods wagons), separation may occur.

Taking the sample just from the surface therefore leads to the wrong result. Representative samples from all layers of the silo can be taken with the SiloDrill.

The SiloDrill fits through all standard openings with its diameter of 90 mm (3.54 in.).



Material	Total length cm (in.)	Spiral screw dia. mm (in.)	Chamber dia. mm (in.)	Volume ml (oz.)	Item no.
Aluminum	150 (59.06)	90 (3.54)	40 (1.57)	400 (13.53)	5301-1000 N4
AISI 304	150 (59.06)	90 (3.54)	40 (1.57)	400 (13.53)	5301-2000 N4



SiloPicker

The SiloPicker is ideal for taking bulk goods samples from silos. Using the extension rods, samples can be taken from depths down to 3.5 m (11.48 ft.). The insertion depth depends on the density of the bulk goods.

▶ Stainless steel AISI 316L (1.4404)



Chamber volume ml (oz.)	Outer dia. mm (in.)	Length cm (in.)	Item no.
100 (3.38)	50 (1.97)	130 (51.18)	5340-1100 N4



Ice borer | for deep-frozen and semi-solid substances

Sampler for deep-frozen products in the food > Stainless steel AISI 316L (1.4404) / AISI and luxury items sector and for semi-solid substances. The very robust hollow drill cuts quickly into the sample material with its sharp

Cross-section sample sawing crown. A sample cylinder (cross-section sample) can be cut out.

- 316 Ti (1.4571)
- Sterilizable



Inner dia. mm (in.)	Outer dia. mm (in.)	Drilling depth mm (in.)	Item no.
1 8 (0.31)	10 (0.39)	100 (3.94)	5323-0190 N4
2 16 (0.63)	21 (0.83)	200 (7.87)	5323-0200 N4
3 50 (1.97)	54 (2.13)	200 (7.87)	5323-2010 N4



Meat sampler for microbiological sampling of carcass surfaces. The samplers for carcass sampling are similar to a cork drill. They can be used in the punch process (a destructive process, also called the cork drill or excision process) to remove a circular sample of tissue of 5 cm² (0.76 in²) from the surface. After the punch sample is cut out, it is removed from the carcass with a scalpel and tweezers.

- For sampling in accordance with DIN EN ISO 17604:2015-12
- Cutting surface 25 mm (0.98 in.) Ø

BeefSteaker



The BeefSteaker meat sampler for manually cutting out a piece of tissue for microbiological sampling of carcass surfaces. The core bit is drilled in by hand.

Drill head



The drill head for carcass sampling has a hexagonal socket for bit holders or chucks, suitable for all commonly used cordless screwdrivers or drilling machines.

Dia. mm (in.)	Item no.	
25 (0.98 in.)	5388-1000	N4

Dia. mm (in.)	Item no.		
25 (0.98 in.)	5388-1010 N4		

Burkle Inc.

155 Keyland Court | Bohemia, NY 11716 Phone +1-631-750-3772 Fax +1-631-750-3773 info@burkle-inc.com | www.burkle-inc.com



Food scoop, blue

Ergonomics were particularly important to us for the design of the sample scoop. The scoop stands horizontally due to the special handle design and the slightly angled front edge and can therefore be put down without the sample material trickling out.

The food scoop is ideal for powders, granulates and pastes.

- PS, blue
- Individually packaged and sterilized
- Production and packaging according to clean room class 7
- Can be used as part of HACCP/IFS foreign object management





Capacity ml (oz.)	Length mm (in.)	Pieces per pack	Item no.
25 (0.85)	141 (5.55)	10	5378-3001 N4
50 (1.69)	160 (6.30)	10	5378-3003 N4
100 (3.38)	205 (8.07)	10	5378-3005 N4
150 (5.07)	216 (8.50)	10	5378-3007 N4
250 (8.45)	232 (9.13)	10	5378-3013 N4
500 (16.91)	280 (11.02)	10	5378-3015 N4
1000 (33.81)	332 (13.07)	10	5378-3019 N4

Measuring spoons, blue

Volumetric and dosing spoons can be used to measure small amounts of powder, granulates and liquids exactly and dose them precisely. For precise measuring, excess medium is stripped off the upper edge of the spoon.

- Individually packaged and sterilized
- Production and packaging according to clean room class 7
- Can be used as part of HACCP/IFS foreign object management



Cap. ml (oz.)	Model	Pcs./ pack	Item no.	
0,5 (0.02)	large package	100	5378-4044	N4
0,5 (0.02)	sterilized, packaged individually	100	5378-4045	N4
1 (0.03)	large package	100	5378-4046	N4
1 (0.03)	sterilized, packaged individually	100	5378-4047	N4
2,5 (0.08)	large package	100	5378-4048	N4
2,5 (0.08)	sterilized, packaged individually	100	5378-4049	N4
5 (0.17)	large package	100	5378-4050	N4
5 (0.17)	sterilized, packaged individually	100	5378-4051	N4
10 (0.34)	large package	100	5378-4052	N4
10 (0.34)	sterilized, packaged individually	100	5378-4053	N4
15 (0.51)	large package	100	5378-4054	N4
15 (0.51)	sterilized, packaged individually	100	5378-4055	N4
25 (0.85)	large package	100	5378-4056	N4
25 (0.85)	sterilized, packaged individually	100	5378-4057	N4
50 (1.69)	large package	100	5378-4058	N4
50 (1.69)	sterilized, packaged individually	100	5378-4059	N4

Spoon for foodstuffs, blue

Available in two sizes, 2.5 ml (0.08 oz.) (similar to a teaspoon) and 10 ml (0.34 oz.) (similar to a tablespoon).

- PS, blue
- Individually packaged and sterilized
- Production and packaging according to clean room class 7
- Can be used as part of HACCP/IFS foreign object management

Length mm (in.)	Capacity ml (oz.)	Pieces per pack	Item no.
127 (5.00)	2.5 (0.08)	100	5378-3011 N4
170 (6.69)	10 (0.34)	100	5378-3012 N4





Single-use tweezers, blue

The disposable forceps are suitable for precise lab work and are ideal for grasping, lifting and moving small objects.

- PS, blue
- Individually packaged and sterilized
- Production and packaging according to clean room class 7
- Can be used as part of HACCP/IFS foreign object management



Shape	Pieces per pack	Item no.	
pointed/angled	100	5378-1050	N4
pointed	100	5378-1051	N4
broad tip	100	5378-1052	N4

Scraper for foodstuffs, blue

The scraper is suitable for scraping off, smoothing or cleaning as well as removing residue.

- PS, blue
- Individually packaged and sterilized
- Production and packaging according to clean room class 7
- Can be used as part of HACCP/IFS foreign object management



Spatula for foodstuffs, blue

As a result of the long, ergonomically-shaped, stable handle and a sharp blade, it is possible to even penetrate directly into containers such as paper or plastic sacks.

- PS, blue
- Individually packaged and sterilized
- Production and packaging according to clean room class 7
- Can be used as part of HACCP/IFS foreign object management



		,				
L x W mm (in.)	Pieces per pack	Item no.	Total length mm (in.)	Insertion depth mm (in.)	Pie	
200 x 80 (7.87 x 3.15)	10	5378-3030 N4	263 (10.35)	150 (5.91)	10	